




Accessories related FAQ: Choose the right knife holder

Knife Holder E	Knife Holder N	Knife Holder NZ	Knife Holder E-TC	Knife Holder GD
High & low profile disposable blades	Steel knives, profile c and d (TC knives)	Tungsten carbide, profile c and d knives (steel knives)	TC 65 disposable blades	Glass knives Diamond knives

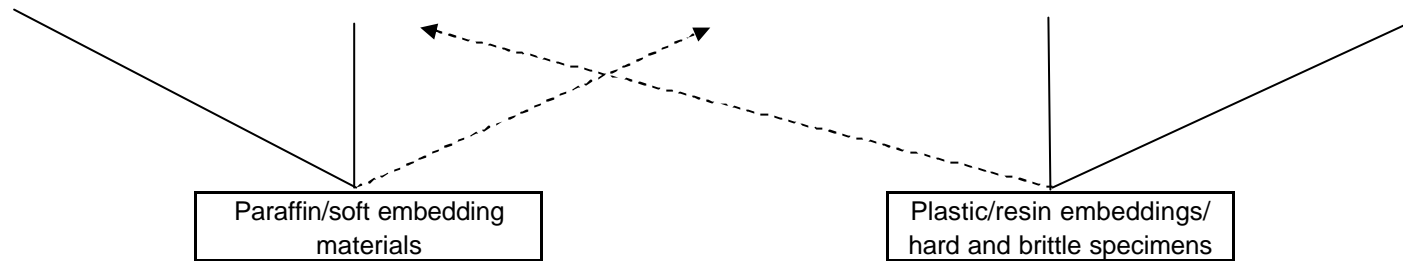

Choose the right specimen clamp

Universal cassette clamp	Standard specimen clamp	Foil clamp	Round specimen clamp
Cassette embedded samples/paraffin	Resin embedded specimens (round/square)	Foils/sheets of plastic	Untrimmed specimens* (mostly industrial specimens)

square round

```

    graph TD
      UCC[Universal cassette clamp] --> CES[Cassette embedded samples/paraffin]
      SSC[Standard specimen clamp] --> RES[Resin embedded specimens (round/square)]
      SSC --> FC[Foil clamp]
      FC --> FSP[Foils/sheets of plastic]
      FC --> RSC[Round specimen clamp]
      RSC --> US[Untrimmed specimens* (mostly industrial specimens)]
  
```

* Trim specimens first to fit them in standard clamp or do a base mounting first to utilize the standard clamp or round specimen clamp.